



Oxygen Radical Absorbance Capacity (ORAC) Assay

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

Tradename: AC Griffonia Lysate Advanced

Code: 16634

CAS #: 999999-99-4

Test Request Form #: 395

Lot #: 27203

Sponsor: Active Concepts, LLC; 107 Technology Drive Lincolnton, NC 28092

Study Director: Erica Segura

Principle Investigator: Meghan Darley

Test Performed:

Oxygen Radical Absorbance Capacity (ORAC)

Introduction

Reactive oxygen species (ROS) are generated by normal cellular processes, environmental stresses, and UV irradiation. ROS are detrimental to cellular structures and functional molecules (i.e DNA, proteins, lipids) as they act as strong oxidizing agents or free radicals. The oxygen radical absorbance capacity (ORAC) assay is a standard method used to assess antioxidant capacity of physiological fluids, foods, beverages, and natural products. The assay quantitatively measures a sample's ability to quench free radicals that have the potential to react with and damage cellular components.

Oxygen Radical Absorbance Capacity (ORAC) assay was conducted to assess the antioxidant capacity of **AC Griffonia Lysate Advanced**.

Assay Principle

This assay is based upon the effect of peroxy radicals generated from the thermal decomposition of 2, 2'-azobis-2-methyl-propanimidamide dihydrochloride (AAPH) on the signal intensity from the fluorescent probe, fluorescein, in the presence of an oxygen radical absorbing substance. The degree of change is indicative of the amount of radical damage and the presence of antioxidants results in an inhibition in the free radical damage to the fluorescein. The antioxidant protection of the sample can be calculated by comparing it to a set of known standards. Trolox®, a water soluble vitamin E analog, with known antioxidant capabilities is used in this ORAC assay as the standard for measuring the antioxidant capacity of unknown substances. ORAC values, expressed in μM of Trolox® equivalents (TE), are calculated using the area under the curves (AUC) of the test product, Trolox®, and the control materials. Trolox equivalency is used as the benchmark for antioxidant capacity of mixtures since it is difficult to measure individual components.

Information contained in this technical literature is believed to be accurate and is offered in good faith for the benefit of the customer. The company, however, cannot assume any liability or risk involved in the use of its chemical products since the conditions of use are beyond our control. Statements concerning the possible use of our products are not intended as recommendations to use our products in the infringement of any patent. We make no warranty of any kind, expressed or implied, other than that the material conforms to the applicable standard specification.



Oxygen Radical Absorbance Capacity (ORAC) Assay

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

Materials

- A. Equipment:** Synergy H1 Microplate reader (BioTek Instruments, Winooski, VT); Gen5 software (BioTek Instruments, Winooski, VT); Pipettes
- B. Buffers:** 75mM Potassium Phosphate (pH 7.4); Deionized H₂O
- C. Reagents:** 2,2'-Azobis(2-methylpropionamidine) dihydrochloride (AAPH) (153mM); 6-Hydroxy-2,5,7,8-tetramethylchromane-2-carboxylic acid (Trolox®); Fluorescein Sodium Salt (4nM)
- D. Preparation:** Pre-heat (37°C) Synergy H1 Microplate reader; Prepare Trolox® standards, sample dilutions, fluorescein solution, and AAPH.
- E. Microtitre Plates:** Corning 96 Well Black Side/Clear Bottom Microplates

Methods

Solutions of **AC Griffonia Lysate Advanced** and Trolox® (positive control) were prepared in 75mM potassium phosphate buffer. Materials were prepared at three different concentrations/dilutions. Trolox® was used as a reference for antioxidant capacity and prepared at a concentrations ranging from 12.5µM to 200µM in 75mM potassium phosphate buffer.

For the ORAC assay, 25µL of test material and Trolox® were combined with 150µL of fluorescein in 75mM potassium phosphate buffer and incubated in the Synergy HT Microplate reader at 37°C for 30 minutes. At the end of the incubation period, 25µL of AAPH were pipetted into each well. Fluorescent measurements were then taken every 2 minutes for approximately 2 hours at an excitation wavelength of 485nm and an emissions wavelength of 520nm.

The AUC and Net AUC values of the standards and samples were determined using Gen5 2.0 Data Reduction Software using the below equations:

$$AUC = 0.5 + \frac{R_2}{R_1} + \frac{R_3}{R_1} + \frac{R_4}{R_1} + \dots + \frac{R_n}{R_1} \rightarrow \text{Where } R \text{ is fluorescence reading}$$

$$Net\ AUC = AUC_{sample} - AUC_{blank}$$

The standard curve was obtained by plotting the Net AUC of different Trolox® concentrations against their concentration. ORAC values of samples were then calculated automatically using the Gen5 software to interpolate the sample's Net AUC values against the Trolox® standard curve. ORAC measurements for the test material were expressed in micro moles Trolox® equivalents (µMTE), where 1 ORAC unit is equal to 1 µMTE.



Results

AC Griffonia Lysate Advanced began exhibiting antioxidant activity at a 0.05% concentration.

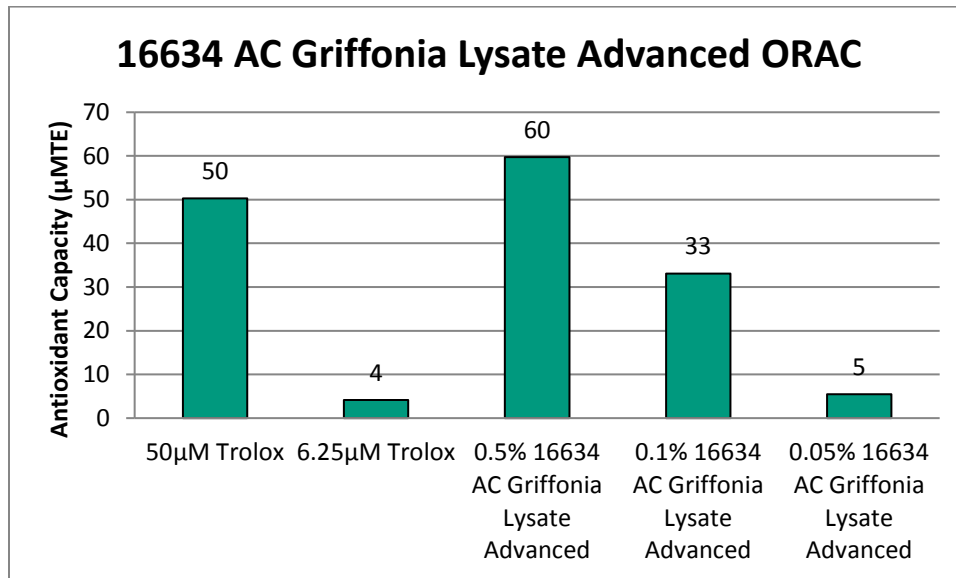


Figure 1: Antioxidant capacities

Discussion

As shown in figure 1, **AC Griffonia Lysate Advanced** exhibited antioxidant activity comparable to 50μM Trolox®. The antioxidant capacity of **AC Griffonia Lysate Advanced** increased as the concentration increased, as a result we can assure that its ability to minimize oxidative stress is dose dependent.

AC Griffonia Lysate Advanced was designed to provide anti-aging, tightening, and lifting properties. With the present study we can confirm that this unique ingredient is not only capable of providing functional benefits but it is also capable of providing potent antioxidant benefits when added to cosmetic applications.