



Phototoxicity Assay Analysis

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

Tradename: PhytoCycle® Orange

Code: 16925

CAS #: 7732-18-5 & 84012-28-2 (or) 8028-48-6 & 1686112-36-6 (or) 68333-16-4

Test Request Form #: 7829

Lot #: N210201E

Sponsor: Active Concepts, LLC; 107 Technology Drive Lincolnton, NC 28092

Study Director: Maureen Danaher

Principle Investigator: Daniel Shill

Test Performed:

In Vitro EpiDerm™ Model (EPI-200-SIT) Phototoxicity

SUMMARY

In vitro phototoxicity irritation studies were conducted to evaluate whether **PhytoCycle® Orange** would induce phototoxic irritation in the EpiDerm™ model assay.

The product was tested according to the manufacturer's protocol. The test article solution was found to be a non-photoirritant at concentrations of 0.5%, 1.5%, 5.0% and 10.0%. Reconstructed human epidermis was incubated in growth media for one hour to allow for tissue equilibration after shipping from MatTek Corporation, Ashland, MA. Test substance was applied to the tissue inserts in four varying concentrations and incubated overnight at 37°C, 5% CO₂, and 95% relative humidity (RH). The following day, the appropriate tissue inserts were irradiated (UVA) for 60 minutes with 1.7 mW/cm² (=6 J/cm²). After substance incubation, irradiation, and washing was completed, the cell viability test was conducted. Cell viability was measured by dehydrogenase conversion of MTT [(3-4,5-dimethyl thiazole 2-yl)], present in the cell mitochondria, into blue formazan salt that was measured after extraction from the tissue. The photoirritation potential of the test chemical was dictated by the reduction in tissue viability of UVA exposed tissues compared to non-UVA exposed tissues.

Under the conditions of this assay, the test article was considered to be non-phototoxic at concentrations of 0.5%, 1.5%, 5.0%, and 10.0%. The negative and positive controls performed as anticipated.

Information contained in this technical literature is believed to be accurate and is offered in good faith for the benefit of the customer. The company, however, cannot assume any liability or risk involved in the use of its chemical products since the conditions of use are beyond our control. Statements concerning the possible use of our products are not intended as recommendations to use our products in the infringement of any patent. We make no warranty of any kind, expressed or implied, other than that the material conforms to the applicable standard specification.



Phototoxicity Assay Analysis

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

I. Introduction

A. Purpose

In vitro dermal phototoxicity study was conducted to evaluate whether a test article would induce photoirritation in the EpiDerm™ model assay. MatTek Corporation's reconstructed human epidermal model is becoming a standard in determining the phototoxicity potential of a test substance. This assay is able to discriminate between photoirritants and non-photoirritants at varying concentrations.

II. Materials

- A. Incubation Conditions:** 37°C at 5% CO₂ and 95% relative humidity
- B. Equipment:** Forma humidified incubator, ESCO biosafety laminar flow hood, Synergy HT Microplate reader; UVA-vis Irradiation Equipment; UVA meter; Pipettes
- C. Media/Buffers:** Dulbecco's Modified Eagle Medium (DMEM) based medium; Dulbecco's Phosphate-Buffered Saline (DPBS); sterile deionized H₂O
- D. Preparation:** Pre-incubate (37°C) tissue inserts in assay medium; Place assay medium and MTT diluent at 4°C, MTT concentrate at -20°C, and record lot numbers of kit components
- E. Tissue Culture Plates:** Falcon flat bottom 96-well, 24-well, and 6-well tissue culture plates
- F. Reagents:** MTT (3-4,5-dimethyl thiazole 2-yl) (1.0mg/mL); Extraction Solution (Isopropanol); Chlorpromazine; Triton X-100 (1%)
- G. Other:** Wash bottle; sterile disposable pipette tips; Parafilm; forceps

III. Test Assay

A. Test System

The reconstructed human epidermal model, EpiDerm™ consists of normal human-derived epidermal keratinocytes which have been cultured to form a multilayer, highly differentiated model of the human epidermis. This model consists of organized basal, spinous, and granular layers, and contains a multilayer stratum corneum containing intercellular lamellar lipid layers. The EpiDerm™ tissues are cultured on specially prepared cell culture inserts.

B. Negative Control

Sterile deionized water is used as the negative controls for the EpiDerm™ Phototoxicity assay.

C. Positive Control

Concentrations of chlorpromazine, ranging from 0.001% to 0.1%, were used as positive controls for the EpiDerm™ Phototoxicity assay.

D. Data Interpretation Procedure

A photoirritant is predicted if the mean relative tissue viability of the 2 tissues exposed to the test substance and 60 minutes of 6 J/cm² is reduced by 20% compared to the non-irradiated control tissues.

Information contained in this technical literature is believed to be accurate and is offered in good faith for the benefit of the customer. The company, however, cannot assume any liability or risk involved in the use of its chemical products since the conditions of use are beyond our control. Statements concerning the possible use of our products are not intended as recommendations to use our products in the infringement of any patent. We make no warranty of any kind, expressed or implied, other than that the material conforms to the applicable standard specification.



Phototoxicity Assay Analysis

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

IV. Method

A. Tissue Conditioning

Upon MatTek kit arrival at Active Concepts, LLC the tissue inserts are removed from their shipping medium and transferred into fresh media and tissue culture plates and incubated at 37°C at 5% CO₂ and 95% relative humidity for 60 minutes. After those 60 minutes the inserts are transferred into fresh media and tissue culture plates and tissue insert dosing begins.

B. Test Substance Exposure

50µL of the diluted test substance in their respective concentrations are applied to 2 tissue inserts and allowed to incubate for overnight in a humidified incubator (37°C, 5% CO₂, 95% RH).

C. Tissue Irradiation

Tissue inserts in their 6-well plates are UVA-irradiated for 60 minutes with 6 J/cm² at room temperature. The non-irradiated tissue inserts are incubated at room temperature in the dark.

D. Tissue Washing and Post Incubation

After UVA-irradiation and dark incubation is complete the tissue inserts are washed using sterile DPBS and transferred to fresh 6-well plates and media for overnight incubation at 37°C, 5% CO₂, 95% RH.

E. MTT Assay

Tissue inserts are transferred into 300µL MTT media in pre-filled plates and incubated for 3 hours at 37°C, 5% CO₂, and 95% RH. Inserts are then removed from the MTT medium and placed in 2mL of the extraction solution. The plate is sealed and incubated at room temperature in the dark for 24 hours. After extraction is complete the tissue inserts are pierced with forceps and 2 x 200µL aliquots of the blue formazan solution is transferred into a 96 well plate for Optical Density reading. The spectrophotometer reads the 96-well plate using a wavelength of 570 nm.

V. Acceptance Criterion

A. Negative Control

The results of this assay are acceptable if the mean negative control Optical Density (OD₅₇₀) is ≥ 0.8.

B. Positive Control

The assay meets the acceptance criterion if a dose dependent reduction in cell viability in the UVA-irradiated tissues is between 0.00316% and 0.0316%.

C. Standard Deviation

Since the phototoxicity potential is predicted from the mean viability of 2 tissues for the EpiDerm™ Phototoxicity Protocol, the variability of the replicates should not exceed 30%.

VI. Results

A. Tissue Characteristics

The tissue inserts included in the MatTek EpiDerm™ assay kit were in good condition, intact, and viable.

B. Tissue Viability Assay

The results are summarized in Figure 1. Cell viability is calculated for each tissue as a percentage of the corresponding vehicle control either irradiated or non-irradiated. Tissue viability was not reduced by 20% in the presence of the test substance and UVA-irradiation at concentrations of 0.5%, 1.5%, and 5.0%. The negative control mean exhibited acceptable relative tissue viability while the positive control exhibited dose dependent loss of tissue viability and cell death.

C. Test Validity

The data obtained from this study met criteria for a valid assay. The negative and positive controls performed as anticipated.

VII. Conclusion

Phototoxicity (photoirritation) is defined as an acute toxic response that is elicited after exposure of the skin to certain chemicals and subsequent exposure to light. Under the conditions of this assay, the test article substance was considered to be non-phototoxic at concentrations of 0.5%, 1.5%, 5.0%, and 10.0%. The negative and positive controls performed as anticipated.

There is a slight decrease in viability at the 10% concentration but viability does not decrease more than the acceptable 20%. We can safely say that **PhytoCycle® Orange** is not a photoirritant when used at the suggested use levels of 1.0% -10.0%.

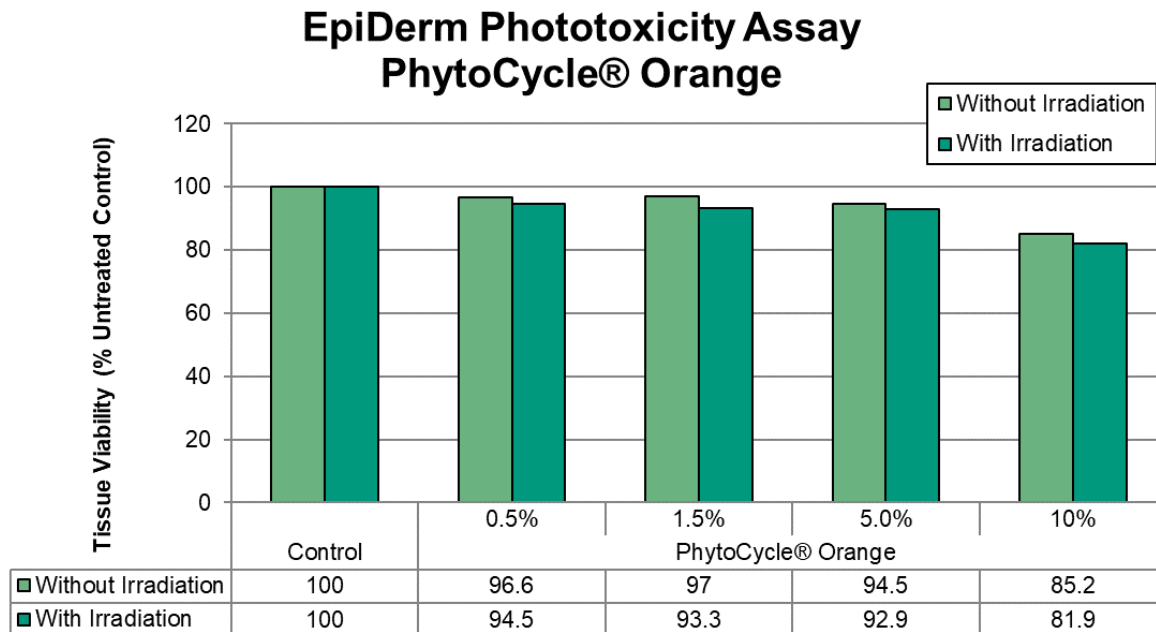


Figure 1: EpiDerm Phototoxicity Graph

Information Contained in this technical literature is believed to be accurate and is offered in good faith for the benefit of the customer. The company, however, cannot assume any liability or risk involved in the use of its chemical products since the conditions of use are beyond our control. Statements concerning the possible use of our products are not intended as recommendations to use our products in the infringement of any patent. We make no warranty of any kind, expressed or implied, other than that the material conforms to the applicable standard specification.