



Cellular Viability Assay Analysis

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

Tradename: ACB Bio-Chelate 5

Code: 20339

CAS #: 7732-18-5 & 8013-01-2 & 8013-01-2 & 8013-01-2 & 8013-01-2 & 8013-01-2

Test Request Form #: 253

Sponsor: *Active Concepts, LLC; 107 Technology Drive Lincolnton, NC 28092*

Study Director: *Erica Segura*

Principle Investigator: *Meghan Darley*

Test Performed:

Cellular Viability Assay

Introduction

The cellular viability assay is useful for quantitatively measuring cell-mediated cytotoxicity, cell proliferation and mitochondrial metabolic activity. Increased metabolism in a cell indicates ample cellular respiration and adenosine triphosphate (ATP) production. ATP is the molecular energy of cells and is required in basic cell function and signal transduction. A decrease in ATP levels indicates cytotoxicity and decreased cell function while an increase in ATP levels indicates healthy cells.

The cellular viability assay was conducted to assess the ability of **ACB Bio-Chelate 5** to increase cellular metabolic activity in cultured dermal fibroblasts.

Assay Principle

The assay utilizes a nonfluorescent dye, resazurin, which is converted to a fluorescent dye, resorufin, in response to chemical reduction of growth medium from cell growth and by respiring mitochondria. Healthy cells that are in a proliferative state will be able to easily convert resazurin into resorufin without harming the cells. This method is a more sensitive assay than other commonly used mitochondrial reductase dyes such as MTT. An increase in the signal generated by resazurin-conversion is indicative of a proliferative cellular state.

<p>This information is presented in good faith but is not warranted as to accuracy of results. Also, freedom from patent infringement is not implied. This information is offered solely for your investigation, verification, and consideration.</p>



Cellular Viability Assay Analysis

info@activeconceptsllc.com • Phone: +1-704-276-7100 • Fax: +1-704-276-7101

Materials

- | | |
|----------------------------------|---|
| A. Kit: | PrestoBlue™ Cell Viability Reagent (Invitrogen, A13261) |
| B. Incubation Conditions: | 37°C at 5% CO ₂ and 95% relative humidity (RH) |
| C. Equipment: | Forma humidified incubator; ESCO biosafety laminar flow hood; Light microscope; Pipettes |
| D. Cell Line: | Normal Human Dermal Fibroblasts (NHDF) (Lonza; CC-2511) |
| E. Media/Buffers: | Dulbecco's Modified Eagle Medium (DMEM); Penicillin-Streptomycin (50U-50mg/mL); Fetal Bovine Serum (FBS); Phosphate Buffered Saline (PBS) |
| F. Culture Plate: | Falcon flat bottom 96-well tissue culture treated plates |
| G. Reagents: | PrestoBlue™ reagent (10X) |
| H. Other: | Sterile disposable pipette tips |

Methods

Human dermal fibroblasts were seeded into 96-well tissue culture plates and allowed to grow to confluency in complete DMEM. A 10-fold serial dilution was performed resulting in **ACB Bio-Chelate 5** concentrations on 1%, 0.1%, and 0.01% in complete DMEM and incubated with fibroblasts for 24 hours.

Ten microliters of viability reagent was added to 90µL of cell culture media in culture wells.

Results

The data obtained from this study met criteria for a valid assay and the controls performed as anticipated.

ACB Bio-Chelate 5 at higher concentrations is able to increase cellular metabolism compared to the control.

Cellular metabolism results are expressed as a percentage of the control.

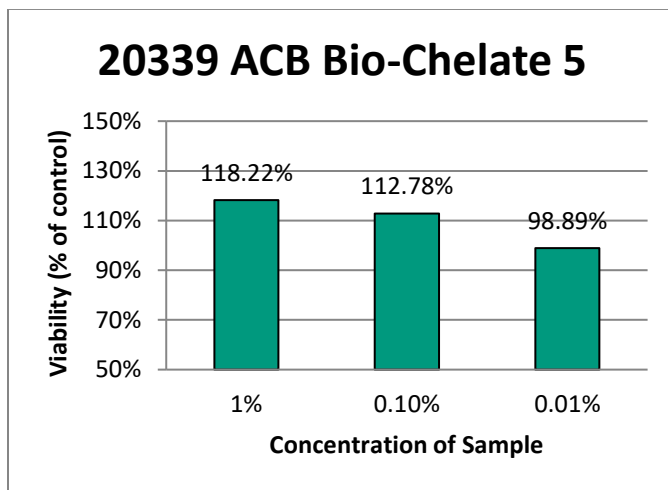


Figure 1: Cellular Metabolism of **ACB Bio-Chelate 5**-treated fibroblasts expressed in terms of percent of control.

Discussion

As shown in figure 1, **ACB Bio-Chelate 5** exhibited positive results by increasing cell metabolism. The increase in fluorescent signal indicates an increase in cellular metabolism and viability post **ACB Bio-Chelate 5** treatment. For these reasons, we can assume **ACB Bio-Chelate 5** is suitable for cosmetic applications designed to increase cell viability and metabolism.